


Rocket® Chelsea & Westminster CVS Set



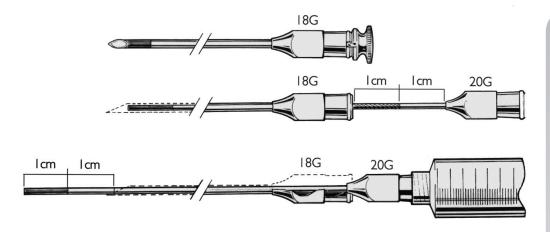

- ◆ THREE PART NEEDLE SET comprising an outer needle with long bevel complete with stylet and an inner cannula needle with blunt distal end to permit secondary sampling. The set is complete with 1cm markings to indicate relative positions of the components
- ♦ HIGH RESOLUTION ECHO BEVEL DESIGN to give smoother, easier tissue penetration with pinpoint ultrasound definition of the needle tip in either amniocentesis and CVS.
- ♦ NEEDLE MARKINGS and depth marker help ensure correct positioning and minimise the risk of over insertion.
- ♦ SECURE LUER LOCK FITTINGS ensure the needle, extension set and syringes remain correctly and securely connected throughout
- ◆ CE MARKED, PURPOSE DESIGNED DEVICE, designed specifically for amniocentesis and transabdominal chorionic villus sampling which avoids the risks of inappropriate use of non-specific generic devices

References:

- 1. Lilford R.J. Chorion villus biopsy. Arch. Dis. Child. 1985 60:897-899
- 2. Maxwell D, Lilford R.J et al. Transabdominal chorionic villus sampling. Lancet Jan 18 1986: 123-126

Description	Order Code
Rocket Chelsea & Westminster pattern CVS Needle Set 18G x 19cm with 20G x 20cm plain ended tube Sterile for single use, packed in carton of 10 sets	R57032

Rocket Medical offer the wellestablished Chelsea & Westminster pattern CVS Needle which includes single and twin sampling functions in the same cost effective needle set.

Indications: For the aspiration of fluid from the amniotic sac and transabdominal chorionic villus sampling for the detection of chromosomal abnormalities in the fetus.

